Tutorial mask-based parameter extraction from images in native space

Scenario:

A mask was generated for each animal in native space. The mask contains two ROIs: ROI-1 with ID-1 (value is 1) and ROI-2 with ID-2 (value is 2).

For a specific image in native space, we want to extract basic parameters such as mean/min/max/number of voxels and volume for each ROI of the mask.

Example: t2.nii-image overlayed with the mask 'EMS.nii', both images in Native space. The mask contains two ROIs.

Content

- 1) Prerequisite
- 2) Steps to perform the parameter extraction
- 2.1) Readout-file selection
- 2.2) Mask-file selection
- 2.3) Select a filename-prefix for resulting excel-file
- 3) The output: Excel-file
- 4) Batch


1) Prerequisite

For all animals the mask file (here: 'EMS.nii') and the readout-file for parameter-extraction (here: 'rc rCBV.nii') exist. Mask-file and readout-file should meet the following criteria:

- share same space (here: native space)
- In register to each other
- Ideally same voxel-size and identical header-matrix

2) Steps to perform the parameter extraction


-Select all animals from the left animal list box

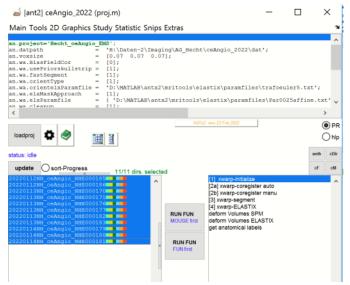

-go to MENU: Statistic/obtain parameter from masks




-this opens the parameter-GUI

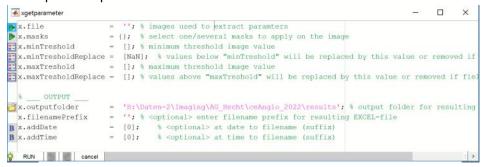

Now do the following steps:


2.1) Readout-file selection

- select green icon left to the **x.file**-parameter
- from the file-section window select "rc_rCBV.nii" → hit [OK].

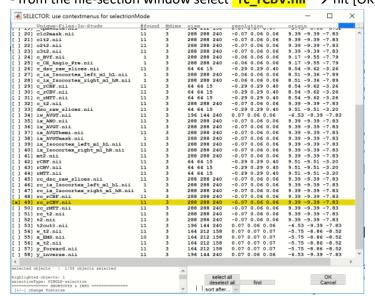


2.2) Mask-file selection

- select green icon left to the x.masks-parameter ...this opens a dialog...

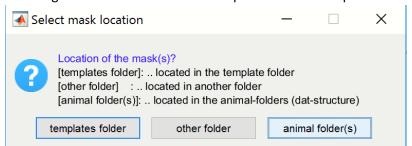

-select [animal folder(s)] from the question-dialog

Reason: In this example each animal has it's own mask (native space), stored in the respective animal-folder (within the "dat"-folder).

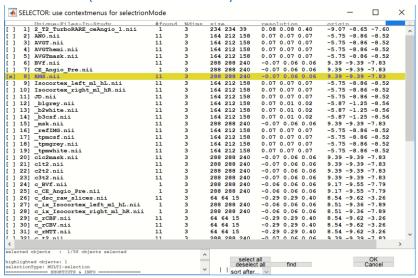

- from the file-section window select the "EMS.nii"-file → hit [OK].

Remember: the "EMS.nii"-file is a mask-file containing 2 ROIs.

Important-1: selection of multiple masks possible here

You can also select more than one mask-image in this step. The numeric-Ids (values) in each mask-file do not interfere with the numeric ids of other selected masks during read-out and creation of the excel-file.

Example: Theoretical scenario: You have 2 masks:

- -Mask1: "mask-1.nii" with 4 ROIs and respective ids: 1,2,3,4 and
- -Mask2: "mask-2.nii" with 3 ROIs and respective ids: 1,2,3
- → You can select the two mask "mask-1.nii" and "mask-2.nii" in one step here
- → The resulting excel-file will contain the mask-file name and the respective ID of the two masks.

Example in the output-file:

- mask-1__ID_1 .. with parameter extraction here mask-1__ID_2 .. with parameter extraction here
- mask-1__ID_3 .. with parameter extraction here
- mask-1__ID_4 .. with parameter extraction here mask-2 ID 1 .. with parameter extraction here
- mask-2__ID_2 .. with parameter extraction here
- mask-2__ID_3 .. with parameter extraction here

Important-2: A text-file with specific labels can be used in combination with the mask

If a text-file with IDs and label is provided, the labels will be used in the resulting excel-file To following is mandatory when using a text-file:

- text-file must have the same name as the mask-file
 - → example: for "EMS.nii" the text-file would by "EMS.txt"
- text-file and mask-file must be in the same folder
 - → in this example both files would be in the animal folders
- the text-file contains: the id (value) and the label for each ROI in the mask (row-wise)
- → for example, for this tutorial the text-file would contain:
 - 1 L somat Cortex
 - 2 R_somat_Cortex

2.3) Select a filename-prefix for resulting excel-file

- for the parameter: **filenamePrefix** type an <u>arbitrary</u> string that is used as prefix in the resulting filename. Here, I used 'EMS_rc_rCBV' (here a combination of mask- and readout-filename) as prefix

The parameter window should now look as follows:

```
xgetparameter
                                                                                            ×
x.file
x.masks
            = {'rc_rCBV.nii' };
                                            % images used to extract paramters
x.masks = {'EMS.nii'}; % select one/several masks to apply on the image x.minTreshold = []; % minimum threshold image value
x.minTresholdReplace = [NaN]; % values below "minTreshold" will be replaced by this value or removed if
  x.maxTreshold = []; % maximum threshold image value
x.maxTresholdReplace = []; % values above "maxTreshold" will be replaced by this value or removed if fiel
      OUTPUT
                     = 'H:\Daten-2\Imaging\AG_Hecht\ceAngio_2022\results'; % output folder for resulting
3x.outputfolder
 x.filenamePrefix = 'EMS rc rCBV'; % <optional> enter filename prefix for resulting EXCEL-file
                     = [0]; % <optional> at date to filename (suffix)
B x.addDate
B x.addTime
                     = [0]; % <optional> at time to filename (suffix)
   RUN 5 cancel
```

Optional Settings

1) <u>minTreshold</u>: set a min-threshold for the readout-image ..., in this case: values below this values will by replace by the value defined by minTresholdReplace (here NAN-values).

The latter can be changed via "minTresholdReplace"-parameter.

- → example: minTreshold = 3 and minTresholdReplace=3: In this case values below 3 will be replaced by the value 3
- 2) maxTreshold/TresholdReplace: ...same logic ...
- 3) outputfolder: change output-folder ...default: "results"-folder in the study-folder
- 4) addDate & addTime: add date & time to the filename of the new excel-file
- -finally hit [RUN]

3. The output: Excel-file


- -When the parameter " $\underbrace{\text{outputfolder}}$ " is not changed the resulting excel-file will be stored in the
- "results"-folder of the current study (aka project).
- -Here the name of the **Excel-file** is "**res_EMS_rc_rCBV_maskParameter.xlsx**" ("res"+filenamePrefix + "_maskParameter.xlsx")

-The Excel-file contains the following sheets: 'info' basic information: date/mask-name/animals/ missing animals (due to missing files) 'frequencyMsk': number of voxels for each ROI (real ROI-size, readout-image independent, true size of the ROI) 'volMsk' : volume of the ROI (real ROI-size, readout-image independent, true size of the ROI) : number of voxels for each ROI (readout-image dependent, 'frequency' NaN/threshold dependent → thus can be smaller than 'frequencyMsk') 'vol' : volume of the ROI (readout-image dependent (Nan-values!)) NaN/threshold dependent → thus can be smaller than 'volMsk') 'mean' : mean intensity value of readout-image within ROI 'std' : standard deviation over intensities of readout-image within ROI 'median' : median intensity value of readout-image within ROI 'integrDens' : product of number voxels and mean intensity of readout-image within ROI 'min' : min intensity value of readout-image within ROI 'max' : max intensity value of readout-image within ROI

- The parameter-sheets have the following structure:

- -ROWS: ...ROI-MASK, i.e., each row represents the parameter readout for a specific ROI
- -COLUMNS:...ANIMALS, i.e. each column contains the data of one animal

Example: mean-parameter (sheet: "mean"): Here we have two ROIs from one mask-image (represented in rows) and a number of animals (columns).

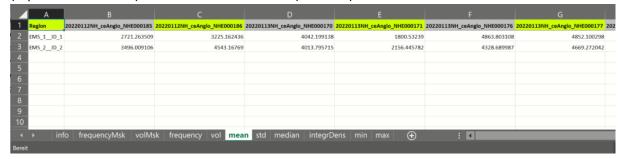

4. Batch

- -type <u>anth</u> to obtain the code... alternatively you could type: **uhelp(anth)** or hit the [anth]-button from the ANT-GUI.
- you can than copy and modify the code

```
FUNCTION:
                                                 [xgetparameter.m]
 z=[].
                                                     'rc_rCBV.nii' };
'EMS.nii' };
                                                                                                                                                                                      images used to extract paramters select one/several masks to apply on the image minimum threshold image value values below "minTreshold" will be replaced by this
                                             = { '
 z.minTreshold
 z.minTresholdReplace = [NaN];
                                                                                                                                                                                    % values below "minTreshold" will be replaced by this
% maximum threshold image value
% values above "maxTreshold" will be replaced by this
% output folder for resulting EXCEL-file
% enter filename prefix for resulting EXCEL-file
% at date to filename (suffix)
% at time to filename (suffix)
 z.maxTreshold
 z.maxTresholdReplace = z.outputfolder = z.filenamePrefix =
                                             = [];
= 'H:\Daten-2\Imaging\AG_Hecht\ceAngio_2022\results';
= 'EMS_rc_rCBV';
 z.addDate
 z.addTime
 xgetparameter(1,z);
₩ <> EM 360/408 #1 ant.m
                                                                        v | 15 | 35 | 70
```

Here is the code to rerun the task: